



Frontmatter

Peripherals

Introduction

Troubleshooting

Basic Physics

Verification

Scanning

Installation

FRUs

Upgrading

Architecture

Advd. Diags

Power System

Service Programs

Configuration

Stress Echo

Diagnostics

Boot Sequence

STB SEQ-010

**acuson**

A Siemens Company


# Sequoia® Ultrasound System

## Service Training Manual



---

## ***Sequoia<sup>®</sup> Ultrasound System Service Training Manual***

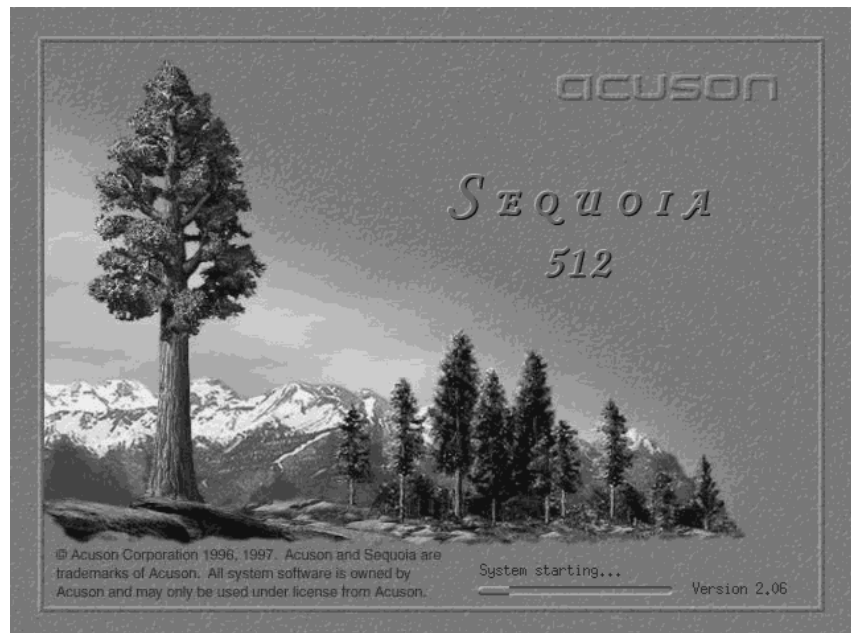

Service Training Department  
Acuson Corporation  
1220 Charleston Road  
Mountain View, CA 94039  
(650) 969-9112  
(800) 4-ACUSON

Manual Assy P/N 58323, Rev. 1  
Language: English  
CE

## COPYRIGHT

Copyright © 2000 by Acuson Corporation. All rights reserved. No part of this publication may be reproduced, transmitted, transcribed, stored in retrieval systems, or translated into any language or computer language, in any form or by any means, electronic, mechanical, magnetic, optical, chemical, manual, or otherwise, without the prior written permission of Acuson Corporation.

Acuson Corporation reserves the right to change its products and services at any time. In addition, this manual is subject to change without notice. Acuson Corporation welcomes customer input on corrections and suggestions for improvements to this manual.

Although Acuson Corporation has attempted to ensure accuracy throughout this manual, Acuson Corporation assumes no liability for any errors or omissions, nor for any damages resulting from the application or use of this information.

## TRADEMARKS

Acuson, AEGIS, CasePro, DELTA, MultiHertz, RES, Sequoia, SpaceTime, The Value of Vision, Vector, 128XP, and XP are registered trademarks of Acuson Corporation. Aspen, Cadence, Convergent, DIMAQ, DTI, EF, FreeStyle, microCase, Multizone, Native, NewView, Perspective, ProtoCALL, QuantX, Quik-Clip, Solo, SST, and WebPro are trademarks of Acuson Corporation.

Cidex, Cidex Plus, and Cidex 7 are registered trademarks of Surgikos, Inc. K-Y Lubricating Gel is a trademark of Johnson & Johnson Products, Inc. IBM is a registered trademark of International Business Machines Corporation. Metricide is a trademark of Metrex Research Corporation. Omnicide is a trademark of Cottrell, Ltd. Panasonic is a trademark of Matsushita Electric Industrial Co. Ltd. Polaroid is a registered trademark of Polaroid Corporation. 3M is a registered trademark of Minnesota Mining and Manufacturing. Apple, AppleTalk, EtherTalk, LocalTalk, Macintosh, MultiFinder, and LaserWriter are trademarks of Apple Computer, Inc. DOT is a registered trademark of Digital Optical Technologies. Pinnacle Micro Inc., Kodak, and Ektascan are registered trademarks of Eastman Kodak Company. PostScript is a registered trademark of Adobe Systems, Incorporated. Ricoh is a trademark of Ricoh Company, Ltd. Sony is a registered trademark of Sony Corporation of America. Verbatim is a registered trademark of Verbatim Corporation. Tosoh is a trademark of Tosoh Corporation. Multi-Imager is a trademark of International Imaging Electronics.

All other product names are trademarks of their respective companies.

**CAUTION!** The Food and Drug Administration of the United States Government restricts use of this device to or on the order of a physician.

## **LICENSE AGREEMENT**

All computer programs copyright 1990-2001 by Acuson Corporation OR ITS SUPPLIERS. Such programs are licensed under the following software license agreement:

Acuson or its suppliers retain(s) ownership of and title to any computer program supplied with the equipment and to the trade secrets embodied in such computer programs. Subject to the Buyer's acceptance and fulfillment of the obligations in this paragraph, Acuson grants the Buyer a personal, non-transferable, perpetual, non-exclusive license to use any computer program supplied with the Equipment that is necessary to operate the Equipment solely on the medium in which such program is delivered for the purpose of operating the equipment in accordance with the instructions set forth in the operator's manuals supplied with the Equipment and for no other purpose whatsoever. Buyer may not reverse-assemble, reverse-compile or otherwise reverse-engineer such computer programs nor may Buyer make a copy of such program or apply any techniques to derive the trade secrets embodied therein. In the event of a failure by Buyer to comply with the terms of this license, the license granted by this paragraph shall terminate. Further, because unauthorized use of such computer programs will leave Acuson without an adequate remedy at law, Buyer agrees that injunctive or other equitable relief will be appropriate to restrain such use, threatened or actual. Buyer further agrees that (i) any of Acuson's suppliers of software is a direct and intended beneficiary of this end-user sublicense and may enforce it directly against Buyer with respect to software supplied by such supplier, and (ii) NO SUPPLIER OF ACUSON SHALL BE LIABLE TO BUYER FOR ANY GENERAL, SPECIAL, DIRECT, INDIRECT, CONSEQUENTIAL, INCIDENTAL OR OTHER DAMAGES ARISING OUT OF THE SUBLICENSE OF THE COMPUTER PROGRAMS SUPPLIED WITH THE EQUIPMENT.

## **REVISION HISTORY**

<b>QRC</b>	<b>P/N-REVISION</b>	<b>INITIATOR</b>	<b>APPROVAL</b>	<b>DATE</b>	<b>CHANGE</b>
			S. Williams	July 1999	Incorporate reviewer comments
A3210	P/N 59191 Rev. 1	J. Madarasz	S. Williams	Dec. 2000	Initial Release.

---

# CONTENTS

<b>MODULE 1</b>	<b>Introducing Sequoia® Ultrasound Technology</b>	
	Sequoia Ultrasound Technology	1-3
	Coherent Imageformer	1-4
	DIMAQ™ Integrated Ultrasound Workstation	1-5
	Doppler Technologies	1-6
	Transducer Technology	1-7
	Worksheet: System Technology	1-8
<hr/>		
<b>MODULE 2</b>	<b>Basic Ultrasound Physics</b>	
	Introducing Ultrasound	2-3
	Ultrasound Applications	2-4
	Ultrasound Imaging	2-5
	Transducers	2-6
	Ultrasound Beam Focusing	2-11
	Ultrasound Wave Reception	2-13
	2-D Artifacts and Enhancements	2-15
	2-D Image Acquisition	2-19
	Scan Conversion	2-20
	Components of Image Quality	2-22
	Other 2-D Features	2-25
	Doppler Principles	2-26
	Spectral Doppler	2-27
	Understanding Blood Flow	2-30
	Artifacts of Spectral Doppler	2-31
	Color Doppler	2-32
	Color Doppler Modes	2-34
	M-Mode	2-35
	Ultrasound Module Questions	2-37

---

**MODULE 3****Scanning**

<b>Module Overview</b>	<b>3-3</b>
<b>Beginning A Scanning Session</b>	<b>3-5</b>
<b>Carotid Artery 2-D Imaging</b>	<b>3-7</b>
<b>Spectral Doppler Imaging</b>	<b>3-9</b>
<b>Color Doppler Imaging</b>	<b>3-12</b>
<b>Liver &amp; Kidney Imaging</b>	<b>3-15</b>
<b>Cardiology/Heart Imaging</b>	<b>3-17</b>
<b>Image Enhancements</b>	<b>3-22</b>
<b>Data Field Values</b>	<b>3-28</b>
<b>Scanning a Phantom</b>	<b>3-30</b>
<b>Scanning Reference Information</b>	<b>3-31</b>
<b>Using AEGIS®</b>	<b>3-32</b>
<b>Scanning Questions</b>	<b>3-34</b>
<b>128XP® User Interface</b>	<b>3-35</b>
<b>Aspen™ 1.0 to 4.0 System User Interface</b>	<b>3-36</b>
<b>Aspen™ and Sequoia Systems New User Interface</b>	<b>3-37</b>

---

**MODULE 4****Field Replaceable Units (FRUs)**

<b>Overview</b>	<b>4-3</b>
<b>Accessing and Removing FRUs</b>	<b>4-4</b>
<b>Brief Board Descriptions</b>	<b>4-8</b>
<b>System Monitor Assembly</b>	<b>4-11</b>
<b>User Interface FRUs</b>	<b>4-13</b>
<b>Covers and Card Cage Assembly</b>	<b>4-19</b>
<b>AUX and FIZ Modules</b>	<b>4-26</b>
<b>Power Subsystem</b>	<b>4-28</b>
<b>SCSI Devices</b>	<b>4-32</b>
<b>Worksheet: FRU Removal and Re-install</b>	<b>4-34</b>

---

**MODULE 5**

<b>System Architecture</b>	
Overview	5-4
System Architecture	5-5
Coherent Imageformer	5-7
Theory of Operation	5-8
Transmitter Board	5-11
Multiplexer Board	5-12
Receiver Board	5-13
Beamformer Board	5-15
Controller Board	5-16
DIMAQ Integrated Ultrasound Workstation	5-17
Color and Spectral Doppler Board	5-22
BDM Board	5-23
Reconstruction and Display Processor Board	5-24
Input/Output Video Board	5-25
Input/ Output Expansion Board	5-26
Peripheral Interface Controller Board	5-27
Physio Interface Module	5-28
Front Panel Processor Board	5-29
2-D/ M-Mode Signal Flow	5-30
Solo™ Spectral Doppler Signal Flow	5-33
Color Doppler Signal Flow	5-35
ECG/Physio Signal Flow	5-37
DIMAQ System Store and Review	5-38
VCR Playback	5-40
Worksheet: System Architecture	5-41

---

**MODULE 6**

<b>Power Distribution</b>	
Overview	6-3
Power Supply Architecture	6-4
Theory of Operation	6-6
Power Subassemblies	6-10
Accessory Isolated Power	6-14
System Cabling Diagram	6-15
Minimum Configuration Requirements	6-16
Power Usage	6-17
Worksheet: Power Supplies	6-18

---

<b>MODULE 7</b>	<b>Configuring the Sequoia System</b>	
	Overview	7-3
	Customizing System Set-up	7-4
	System Configuration	7-13
	Board Configuration	7-14
	Worksheet: Configuration	7-15
<hr/>		
<b>MODULE 8</b>	<b>Diagnostic Tools</b>	
	Overview	8-3
	Learning MSM Tools	8-4
	Expert System Tools	8-10
	Viewing Sequoia System Logs	8-13
	Worksheet: MSM Tools	8-17
<hr/>		
<b>MODULE 9</b>	<b>Installing Peripherals</b>	
	Questions	9-3
<hr/>		
<b>MODULE 10</b>	<b>Sequoia Troubleshooting Tools</b>	
	Overview	10-3
	Recommended Approach	10-4
	Identifying Symptoms	10-5
	Stage 1 - System Down	10-6
	Stage 2 - No Boot	10-8
	Stage 3 - Functionality Problem	10-10
	Stage 4 - Subtle Image Quality Problems	10-11
	Stage 5 - Lock-ups	10-12
	Stage 6 - Intermittent	10-12
	Choosing Appropriate Tools for the Symptom	10-13
	Managing Sequoia System Logs	10-16
	Worksheet: Fault Isolation	10-20
	Troubleshooting Scenarios	10-21
<hr/>		



<b>MODULE 11</b>	<b>Verifying the Sequoia System</b>	
	Overview	11-3
	Mechanical Inspection	11-4
	Electrical Verification	11-5
	Preparing for Performance Verification	11-6
	Hardware Verification	11-8
	Functional Verification	11-9
	AEGIS System Functional Verification	11-15
	Peripheral Device Verification	11-16
	Worksheet: Verification	11-17
<hr/>		
<b>MODULE 12</b>	<b>Installing the Sequoia System</b>	
	Overview	12-3
	Worksheet: Installation	12-4
<hr/>		
<b>MODULE 13</b>	<b>Upgrading the Sequoia System</b>	
	Overview	13-3
	Worksheet: Upgrades	13-4
<hr/>		
<b>MODULE 14</b>	<b>Advanced Diagnostic Tools</b>	
	Overview	14-3
	Connecting the Laptop Hardware	14-4
	Setting Hyper Terminal	14-5
	Setting Reflection Options	14-6
	Using the Sequoia Terminal	14-8
	Running Diagnostics	14-11
	Saving Data	14-12
	Worksheet: Field Laptop	14-13
<hr/>		
<b>MODULE 15</b>	<b>Acuson's Service Programs</b>	
<hr/>		

<b>MODULE 16</b>	<b>Stress Echo Cardiography</b>	
	Overview	16-3
	Defining Stress Echo Cardiography	16-4
	Exercise Stress Program	16-5
	Pharmacological Stress Program	16-6
	How Does It Work?	16-8
	Connecting the ECG from OEMs	16-13
	Typical Stress Echo Problems and Solutions	16-16
	Worksheet: Stress Echo	16-18
<hr/>		
<b>MODULE 17</b>	<b>Interpreting Boot Sequence</b>	
	Overview	17-3
	Sequoia Boot Sequence	17-4
	Sequoia 3.nn Boot Details	17-6
	Boot Messages	17-7
	Boot Sequence Example	17-10
	Worksheet: Boot Sequence	17-14
<hr/>		

## PREFACE

Acuson is pleased to welcome you to the Sequoia® ultrasound systems. The Sequoia systems are special-purpose machines designed for general radiology, obstetrics, gynecology, cardiology, and/or vascular exams.

Acuson systems have an exceptional record of dependability, and Acuson has a customer service network ready to respond to the needs of our customers.

For information on Acuson's service offerings, please consult the Acuson Service Provider or call us at 800-9-ACUSON or 650-969-9112. For local numbers outside the United States, see the Contacting Acuson heading in this section.

### ABOUT THIS MANUAL

The purpose of this manual is to familiarize service personnel with basic operation, maintenance, and troubleshooting of Sequoia systems. Service personnel are expected to be familiar with the use of basic test equipment (DVM, oscilloscope, etc.) and with basic ultrasound theory.

### WARNING!

This manual is intended for use by trained service personnel. There are lethal voltages inside Sequoia systems. Bodily harm, possibly death, and damage to the system may result from untrained individuals opening the system.

### KEY CONVENTIONS

Use this manual when servicing a Sequoia ultrasound system.

This manual uses several conventions to refer to controls on the system. The following list shows the conventions and their descriptions.

Symbol	Description
[LEFT]	Represents a soft key
<b>CALC</b>	Represents a key, knob, switch, or toggle control
<b>CODE + SIZE</b>	Means "Press the <b>CODE</b> key and press <b>SIZE</b> "
◆	Indicates steps to follow to perform a procedure

### NOTE:

A Note indicates especially important information. It is not an alert to a dangerous situation.

### CAUTION!

A Caution alerts the user to the potential for damage to the ultrasound system or other equipment.

**WARNING!** A Warning alerts the user to a potential risk of death or serious injury to the user or patient.

Acuson provides special alphanumeric keys in several languages; however, in this manual, labels for all keys and annotation terms are given in English. To convert the system to another language, contact the Acuson Customer Service representative. For English and international keyboard layouts, see the Keyboard section of the *User Manual*.

## OTHER SEQUOIA SYSTEM MANUALS

In addition to this manual, the following manuals are available for the Sequoia system:

### Acuson Sequoia® Ultrasound System User Manual

The *User Manual* provides information used regularly to perform ultrasound exams. It explains how to operate the standard and optional features of the Sequoia system for radiology, obstetrics, cardiology, gynecology, and/or vascular exams, including the use of special purpose transducers such as endocavity transducers, transesophageal transducers, and interoperative transducers. It explains how to use transducer needle guides to perform biopsies. It also contains setup information, which is used less frequently. The *User Manual* addresses the reader who is familiar with ultrasound techniques; therefore, it does not include sonography training or clinical procedures.

### Acuson Sequoia® Ultrasound System Reference Manual

The *Reference Manual* consists of the *Safety Manual*, the *Transducer Specifications Manual*, and the *Administrator Manual*.

- The *Safety Manual* explains system and transducer safety precautions and maintenance procedures.
  - The *Transducer Specifications Manual* lists technical information about the individual transducers, including power value tables for all Acuson transducers.
  - The *Administrator Manual* describes system and calculation package setup options. In addition, it contains information on setting up and operating peripheral devices with the Sequoia system.
-

## CONTACTING ACUSON

For additional information about the Acuson system, contact the nearest Acuson office.

### U.S.A.

Acuson Corporation, U.S.A.  
Telephone: 1-(650) 969-9112  
1-800-4-ACUSON  
Fax: 1-(650) 964-8331

### ASIA PACIFIC/LATIN AMERICA

Acuson Asia Pacific/Latin America, U.S.A.  
Telephone: 1-(650) 969-9112  
1-800-4-ACUSON  
Fax: 1-(650) 965-9471

### AUSTRALIA

Acuson Pty. Ltd.  
Telephone: (61) 29201-7777  
Fax: (61) 29869-1600

### BELGIUM

Acuson SA/NV  
Telephone: (32) 761-06-96  
Fax: (32) 761-06-90

### CANADA

Acuson Canada Ltd.  
Telephone: 1-(416) 296-1634  
Fax: 1-(416) 296-1259

### DENMARK

Acuson AS  
Telephone: (45) 39179991  
Fax: (45) 40300796

### EUROPE

Acuson EHQ.  
Telephone: (44) 1895-202960  
Fax: (44) 1895-202965

### GREAT BRITAIN

Acuson Ltd., England  
Telephone: (44) 1895-251010  
Fax: (44) 1895-202999

### FINLAND

Acuson OY  
Telephone: (358) 9-27092290  
Fax: (358) 9-27092299

### FRANCE

Acuson S.A.R.L.  
Telephone: (33) 1-6929-8110  
Fax: (33) 1-6446-9840

### GERMANY

Acuson GmbH  
Telephone: (49) 91-31-82650  
Fax: (49) 91-31-826565

### HONG KONG

Acuson Hong Kong Ltd.  
Telephone: (852) 2574-7989  
Fax: (852) 2574-6945

### ITALY

Acuson S.p.A.  
Telephone: (39)02 982-591  
Fax: (39)02 982-41107

### JAPAN

Acuson Nippon, K.K.  
Telephone: (81) 3-3316-3999  
Fax: (81) 3-3316-7002

### LATIN AMERICA

Acuson Latin America, U.S.A.  
Telephone: 1-(305) 596-3148  
Fax: 1-(305) 596-3149

### NORWAY

Acuson AS  
Telephone: (47) 67-58-02-00  
Fax: (47) 67-58-15-30

### RUSSIA

Acuson Int'l Sales Corp.  
Moscow, Russia  
Telephone: (7) 095-956-0018  
or (7) 095-956-0019  
Fax: (7) 502-222-1269

### SINGAPORE

Acuson Singapore Pte Ltd.  
Telephone: (65) 254-8982  
Fax: (65) 254-8952

### SWEDEN

Acuson AB  
Telephone: (46) 8-594-40660  
Fax: (46) 8-594-40675

### SPAIN

Acuson Iberica S.A.  
Telephone: (34) 91-372-17-55  
Fax: (34) 91-372-14-78

---

## COURSE PROCEDURES

### **COURSE ATTENDANCE**

This class runs Monday through Thursday, from 8:00 A.M. to 5:30 P.M. in the Building H Training Room, with a break for lunch at about 11:30 A.M. In addition, students may stay after class or contact the instructor for additional help as needed.

There will be breaks between the sessions; however, please be sure to return to the training room on schedule.

---

### **MODULES**

Each module includes lab instructions, reference material, and evaluation questions. When the module is complete, confer with the instructor to determine whether the objectives have been met before going on to the next module.

---

### **RESOURCES**

All references for the course are in the course modules and handouts. Lab equipment is at each lab station. Please do not remove supplemental resources from the room. For copies, contact an instructor.

Please feel free to ask questions at any time. There will be Post-It notes by the back door. For an issue or question that cannot be resolved immediately, make a note of it on a Post-It note and stick it on the back door.

---

### **PRACTICE**

Feel free to use any classroom equipment to practice the module skills. Ask for help or clarification as needed.

---

## COURSE SCHEDULE

The schedule is approximate and may vary. Short breaks will be taken between sections as appropriate. One hour will be given for lunch some time between 11:30 A.M. and 1:30 P.M.

### MONDAY

TIME	TOPIC
8:30 A.M.	Course introduction & Acuson history
9:30 A.M.	Introduction to Ultrasound Physics
11:00 A.M.	Scanning lecture and lab: carotid artery
2:00 P.M.	General imaging demo and lab
3:00 P.M.	Cardiology demo and lab
4:00 P.M.	Image Enhancements and AEGIS

### TUESDAY

8:00 A.M.	Replacing field replaceable units (FRUs)
10:30 A.M.	Scanner architecture lecture
1:00 P.M.	Verifying the Scanner lab
2:15 P.M.	System troubleshooting lab - Scanner

### WEDNESDAY

8:00 A.M.	System architecture- DIMAQ workstation
9:30 A.M.	Diagnostics lecture and lab
11:00 A.M.	Software architecture lecture
1:15 P.M.	Factory tour
2:00 P.M.	System troubleshooting- DIMAQ workstation

### THURSDAY

8:00 A.M.	Configuration lecture and lab
9:20 A.M.	Power supply lecture and lab
11:00 A.M.	Troubleshooting lecture
11:30 A.M.	Acuson luncheon
1:00 P.M.	Power supply and whole system troubleshooting
4:00 P.M.	Verification lecture and lab
5:00	Course conclusion: certificates and evaluations